



APPLICATION NOTE

Subject : Standard sample analysis of Vitamin A

▶ Analysis condition


Mobile phase : A=Acetonitril, B=MeOH (5:2)

Flow rate : 1ml/min

Detection wavelength : 325nm

Sample injection : 20uL

▶ Chromatogram

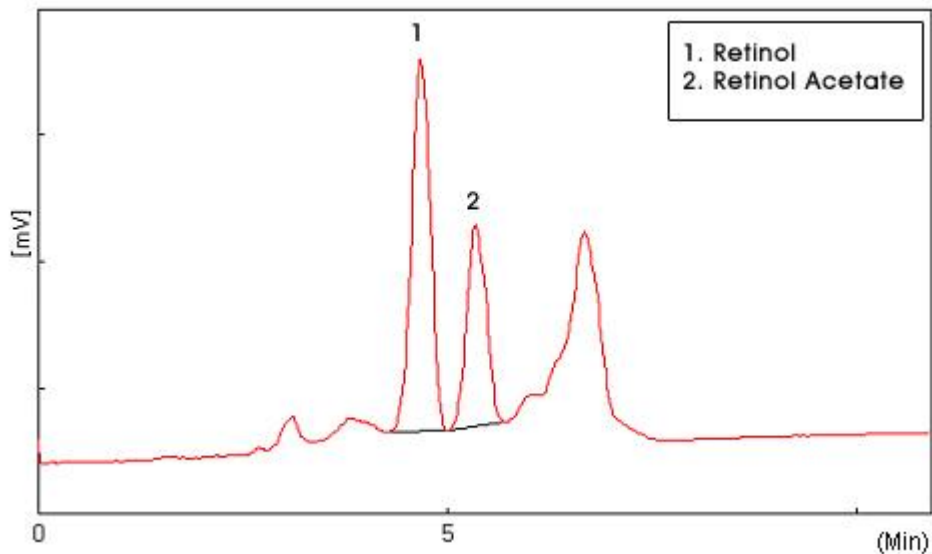

<Fig 1. Analysis of Vitamin A>



영린기기

431-836 경기도 안양시 동안구 호계2동 899-6 Tel : 031-428-8700 Fax : 031-428-8787
Homepage : www.younglin.com E-mail : sales@younglin.com